

ScubaMind - Mental Health

Mental Health Monitoring Application

A privacy-preserving mobile application for behavioral risk awareness, goal tracking, AI-supported guidance, and local-first wellbeing insights.

USER MANUAL



Project Website & Online Resources: <https://scuba-mental-health.netlify.app/>

Version 0.1.0 | Bilkent University — Senior Design Project
Team T2507 | May 2026

Team Members

Merve Güleç • Metin Çalışkan • Muhammed Fatih Başal • Murathan Işık • Yiğit Koşum

Supervisor: Hamdi Dibeklioglu

Table of Contents

1 Introduction	4
1.1 Purpose of This Manual	4
2 Quick Start Guide	4
3 System Requirements	5
4 Installation	5
4.1 Installing from APK	5
4.2 Installing Health Connect (required for wearable data)	6
5 First Launch & Onboarding	6
5.1 Welcome Screen (Step 1)	6
5.2 PHQ-8 Mental Health Questionnaire (Steps 2–10)	6
5.3 Requesting Permissions	8
5.4 Optional: Enabling the Specialized Keyboard	8
5.5 Optional: Samsung Galaxy Watch Integration	9
6 Home Screen	9
6.1 Screen Layout Overview	9
6.2 Mental Risk Level Card	11
6.3 Risk Trend Chart	11
6.4 Today's Highlights	11
6.5 Scuba Recommendations	11
6.6 Logging Your Mood	12
6.7 Pull-to-Refresh	12
7 Goals	12
7.1 Active Goals	15
7.2 Updating Goal Progress	15
7.3 Creating a New Goal	15
7.4 Discover New Goals	16
7.5 Voice Recordings	16
7.6 PHQ-8 Check-in	16
7.7 Daily Voice Check-in	17
7.8 Goal History & Badges	17
7.9 Dive Badges — Achievements	17
8 Chat - Scuby AI Companion	18
8.1 Chat Screen Layout	18
8.2 Starting a Conversation	19
8.3 How Scuby Uses Your Data	19
8.4 Safety Guardrails	20
8.5 Clearing a Conversation	20
9 Profile	20
9.1 User Card	22
9.2 Stats Bar	22
9.3 Daily Check-in	22
9.4 Deep State	22

9.5 Recommendations.....	22
9.6 Achievements.....	23
9.7 Tracking Logic.....	23
9.8 Privacy & Data (Quick Access).....	23
9.9 Global Settings.....	23
9.10 Support.....	23
10 Risk Breakdown.....	23
10.1 Signal Contribution Radar.....	25
10.2 Composite Risk Trend.....	25
10.3 Summary Section.....	25
10.4 Why This Score? — How to Read This.....	25
11 Detailed Signal Pages.....	26
11.1 Activity Page.....	26
11.2 Sleep Page.....	27
11.3 Mood Page.....	28
11.4 Phone Usage Page.....	29
11.5 Mobility Page.....	30
11.6 Mobility Signal (AI Analysis) Page.....	31
11.7 Wearable Signal Page.....	33
12 Privacy & Data.....	34
12.1 Passive Signals.....	34
12.2 Content & Messages.....	35
12.3 Keyboard.....	35
12.4 Research & Analytics.....	35
12.5 Data Controls.....	35
12.6 Complete Data Purge (Delete All Data).....	35
13 Settings.....	36
13.1 General.....	37
13.2 Notifications.....	37
13.3 About.....	37
14 Notifications.....	38
14.1 Notification Categories.....	38
14.2 Managing Notifications.....	38
15 Exporting a Report.....	39
15.1 Generating a PDF Report.....	40
15.2 What the Report Contains.....	41
16 Data Security & Privacy.....	41
16.1 Local-First Architecture.....	41
16.2 AES-256 Encrypted Database.....	42
16.3 Android Application Sandbox.....	42
16.4 Privacy and Security Design Considerations.....	42
17 Troubleshooting & FAQ.....	42
18 Disclaimers & Important Notices.....	43
19 Contact & Feedback.....	44

1 Introduction

ScubaMind is a privacy-preserving mobile mental health monitoring application developed as a Senior Design Project at Bilkent University. It runs on Android smartphones and optional Samsung Galaxy Watch wearables.

ScubaMind continuously and passively collects behavioral and physiological signals — including step counts, sleep patterns, GPS-based mobility, phone usage, voice recordings, and keyboard text — and analyses them using on-device artificial intelligence models. The results are presented to the user as a weekly Behavioral Risk Score along with personalized insights and recommendations.

Important — What ScubaMind Is and Is Not

- ✓ ScubaMind is a decision-support and behavioral awareness tool.
- ✓ It helps you notice patterns in your daily behavior over time.
- ✓ It encourages you to seek professional support when necessary.
- ✗ ScubaMind is NOT a medical or clinical diagnostic tool.
- ✗ Its outputs should never replace a qualified mental-health professional.

All sensitive data is stored and processed exclusively on your device. No raw behavioral data, voice recordings, or keyboard text ever leaves your phone without your explicit consent.

1.1 Purpose of This Manual

This user manual explains how to use ScubaMind from installation to daily usage. It guides users through onboarding, permissions, the Home Dashboard, risk interpretation, goals, notifications, privacy settings, AI chat support, and report generation. It is intended for first-time users, project evaluators, and anyone who wants to understand how ScubaMind works in practice.

2 Quick Start Guide

1. Install ScubaMind on your Android device.
2. Open the application and tap Start on the welcome screen.
3. Complete the PHQ-8 onboarding check-in or skip optional steps when available.
4. Grant the permissions you want to use, such as activity, location, notifications, microphone, usage access, and Health Connect.
5. Reach the Home Dashboard and review your current mental risk level.
6. Tap View breakdown to understand which behavioral signals contribute to your score.
7. Create wellbeing goals and enable reminders if desired.
8. Use Scuby AI Chat for supportive and reflective guidance.
9. Manage data sources from Profile → Privacy & Data.
10. Export a PDF wellbeing report from Privacy & Data → Export your data.

3 System Requirements

Requirement	Minimum / Recommended
Android Version	Android 11 (API level 30) or higher
RAM	At least 3 GB RAM (for on-device ASR inference)
Storage	At least 500 MB of free internal storage
Internet	Required only for the AI Chat (Gemini API). All other features work offline.
Wearable (optional)	Samsung Galaxy Watch (Wear OS) with Health Connect installed
Google Health Connect	Required on the device for sleep & HRV data (available on Play Store)

Offline Use Note

Most ScubaMind features work offline because behavioral data processing and risk analysis are performed locally on the device. Internet access is only required for the Scuby AI Chat feature, which uses the Gemini API.

Battery & Performance

ScubaMind runs background sensing services continuously. To ensure accurate data collection, do not restrict its battery optimization in Android Settings.

The on-device Whisper speech recognition model requires significant CPU resources. Inference may take up to 10 seconds per 30-second audio clip on mid-range devices.

4 Installation

4.1 Installing from APK

ScubaMind is distributed as an Android APK file. It is not yet available on the Google Play Store. Follow the steps below to install it manually.

1. Obtain the ScubaMind APK file from your project supervisor or the team's distribution link.
2. Transfer the APK to your Android device via USB cable, email, or cloud storage.
3. On your device, open Settings → Apps → Special app access → Install unknown apps.
4. Select the app you will use to open the APK (e.g., Files, Chrome) and enable Allow from this source.
5. Locate the APK file using your file manager and tap it.
6. Tap Install when prompted.
7. Once installation completes, tap Open or find ScubaMind in your app drawer.

Security Notice

Enabling 'Install unknown apps' allows installation from sources outside the Play Store. Disable this permission after installation is complete to maintain device security.

Always verify the APK file was received from a trusted source before installing.

4.2 Installing Health Connect (required for wearable data)

If you wish to use wearable physiological data (steps, sleep, HRV from Samsung Galaxy Watch):

8. Open the Google Play Store and search for Health Connect.
9. Install the Health Connect application by Google LLC.
10. Open Health Connect and follow the on-screen setup.
11. When ScubaMind requests Health Connect permissions during onboarding, grant them to enable wearable data integration.
12. After installing ScubaMind, open the app and continue with the onboarding flow. The app will guide you through the initial check-in, permission explanations, optional wearable setup, and the Home Dashboard.

5 First Launch & Onboarding

The first time you open ScubaMind, you will be guided through a 10-step onboarding process. This setup takes approximately 2–12 minutes, depending on the options you choose. You can skip optional steps at any time.

5.1 Welcome Screen (Step 1)

The opening screen displays the ScubaMind branding with an underwater theme. Tap Start to begin onboarding.

What you see on the Welcome Screen

- 'SCUBA Check-in' header with a step counter (Step 1 / 10)
- Estimated time shown in the top-right corner
- Tagline: 'Take a deep breath. Dive into your mind.'
- Description: A 60-second check-in to understand how you've been feeling
- Pill badge: 'Private • Quick • Personalized'
- Note: 'You can skip the optional questions any time.'
- Navigation buttons: Back and Start

5.2 PHQ-8 Mental Health Questionnaire (Steps 2–10)



ScubaMind uses a modified PHQ-8 mental health questionnaire as a baseline check-in tool. Each question asks how often the user has experienced a particular feeling over the past two weeks.

Answer Option	Meaning
Not at all	0 — You have not experienced this
Several days	1 — Experienced on some days
More than half the days	2 — Experienced most days
Nearly every day	3 — Experienced almost every day

Example question shown during onboarding:

Question 4: Feeling tired or having little energy

Tap your answer and then tap Next to proceed. You can go back at any time to change a previous answer.

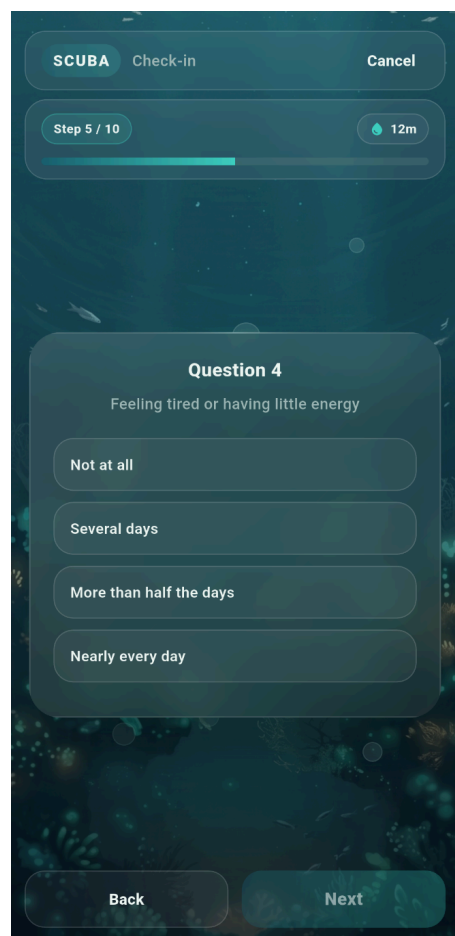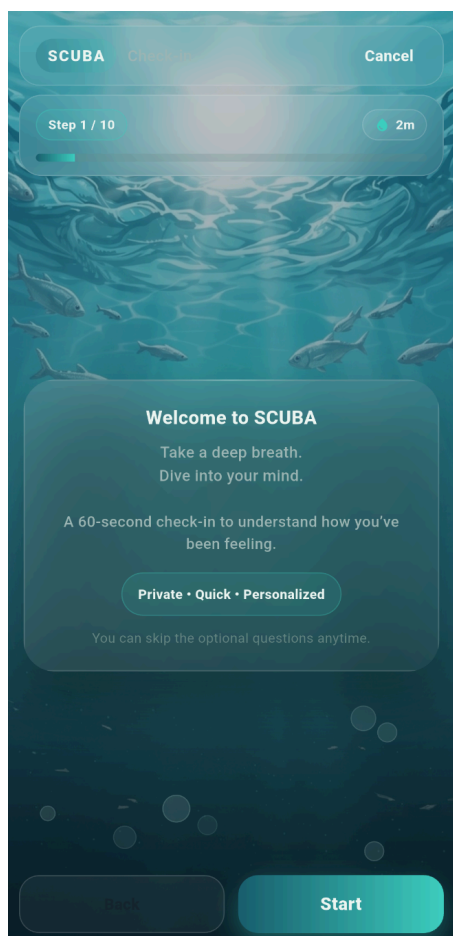

⚠️ PHQ-8 Data Usage

Your PHQ-8 responses are stored only in the encrypted local database on your device. They are used to calibrate the AI model and as a reference label for behavioral analysis. PHQ-8 data is never stored in plain-text SharedPreferences and is never sent to external servers without your explicit opt-in consent for research participation.

5.3 Requesting Permissions

After the questionnaire, ScubaMind requests the permissions needed to collect behavioral data. Each permission is explained before it is requested.

Permission	Purpose
Physical Activity	Step counting via the built-in step counter sensor
Location (precise)	GPS-based mobility tracking — radius of gyration, routine index, places visited
Microphone	Voice journaling and prosody-based behavioral analysis
Usage Access (UsageStats)	Monitoring app usage duration and screen unlock frequency
Notifications	Goal reminders, check-in prompts, and system alerts
Health Connect (optional)	Importing HRV and sleep stage data from Samsung Galaxy Watch

You can grant, deny, or later revoke any permission at any time from Profile → Privacy & Data. Denying permission does not prevent ScubaMind from working — the system gracefully degrades and uses only the available data sources.

Permission Denied: What Happens in ScubaMind		
Denied Permission	What Happens in ScubaMind	
 Physical Activity	Step-based activity insights may be unavailable, but other available signals continue to work.	
 Location	Mobility analysis, routine index, and place diversity cannot be calculated.	
 Microphone	Voice check-ins and voice-based analysis are disabled.	
 Usage Access	Phone usage and screen-time insights cannot be generated.	
 Notifications	Goal reminders and check-in alerts will not appear.	
 Health Connect	Wearable-based sleep, HRV, and activity data will not be imported.	
 Keyboard	Text-based linguistic analysis remains disabled.	

5.4 Optional: Enabling the Specialized Keyboard

ScubaMind includes a custom Input Method Editor (IME) keyboard that can capture anonymized linguistic features from your text across applications. This feature is entirely optional and disabled by default.

To enable it:

13. Go to Android Settings → General Management → Keyboard list and default.
14. Add ScubaMind Keyboard to the list of enabled keyboards.
15. Optionally, set it as your default keyboard.
16. Return to ScubaMind and toggle Keyboard Analysis on in Privacy & Data.

Privacy of Keyboard Data

Raw typed text NEVER leaves your device.

Only anonymized linguistic features (sentiment scores, complexity metrics) are stored locally.

You can disable keyboard analysis at any time from Privacy & Data.

You can select which apps are monitored under Privacy & Data → Keyboard → Keyboard Apps.

5.5 Optional: Samsung Galaxy Watch Integration

If you own a Samsung Galaxy Watch, you can connect it via Health Connect to import Heart Rate Variability (HRV) and sleep stage data into ScubaMind.

17. Ensure Samsung Health and Health Connect are both installed and up to date on your phone.
18. Open Health Connect and confirm that Samsung Health is allowed to write data.
19. During onboarding, when prompted, tap Grant Health Connect Permissions and allow access to sleep, steps, and heart rate data.
20. ScubaMind will automatically sync wearable data through the Health Connect Bridge component.




6 Home Screen

The Home screen is the central dashboard of ScubaMind. It provides an at-a-glance overview of your mental health risk level, daily behavioral highlights, and personalized recommendations. After onboarding, every subsequent app launch takes you directly to this screen.

6.1 Screen Layout Overview

UI Element	Description
App Bar	Displays 'ScubaMind', today's date, and a notification bell icon
Greeting Banner	Shows time-appropriate greeting ('Good afternoon') with animated mascot
Tracking Status	Shows 'Calm tracking on — Keep going, you're building consistency.'
Log Mood Button	Quick shortcut to log your current mood
Mental Risk Level Card	Gauge showing your current risk score (0–100) and qualitative label
Risk Trend (7 days)	Line chart showing your risk score over the past 7 days

UI Element	Description
Today's Highlights	Four insight cards (risk vs yesterday, mood trend, activity status, weekly avg)
Scuba Recommendations	2–3 contextual, AI-generated wellness suggestions
Bottom Navigation Bar	Home / Goals / Chat / Profile

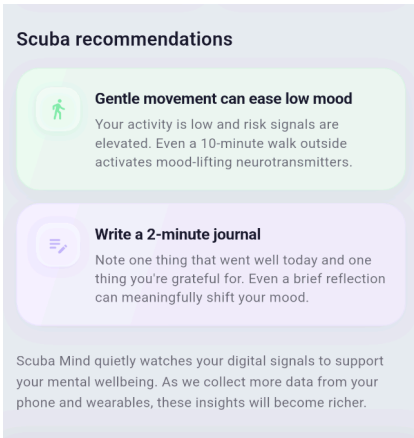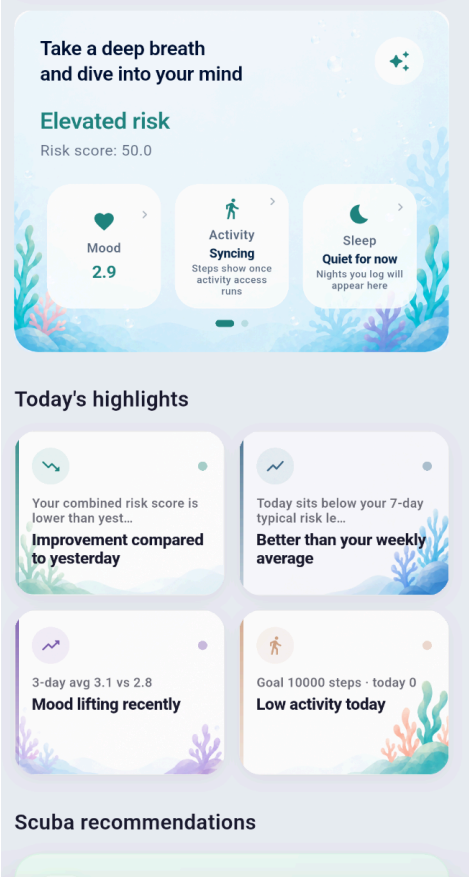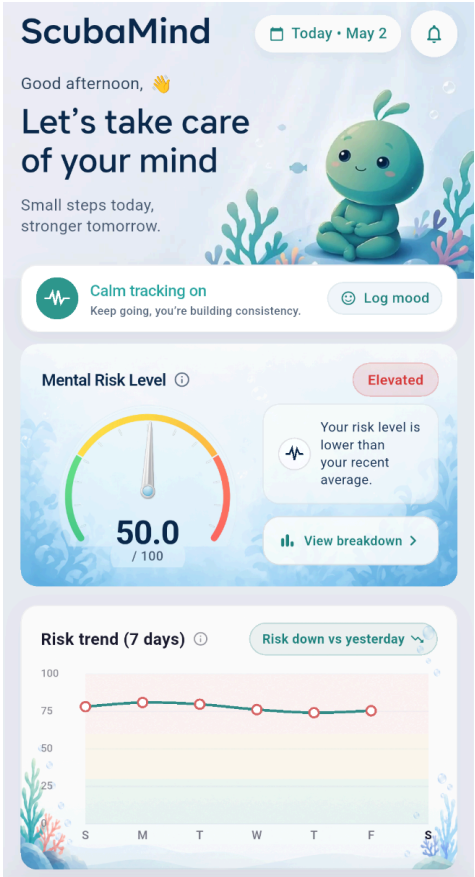

6.2 Mental Risk Level Card

This is the most prominent element on the Home screen. It displays:

- A semi-circular gauge ranging from 0 (low risk) to 100 (high risk), colored green → yellow → red
- Your current composite risk score is a large number (e.g., 50.0 / 100)
- A qualitative risk label: Low, Moderate, Elevated, or High
- A brief contextual message (e.g., 'Your risk level is lower than your recent average.')
- A 'View breakdown >' link to the detailed Risk Breakdown page

Understanding Your Risk Score

The risk score is computed weekly using a rolling 7-day window of multimodal behavioral data. It combines: GPS mobility patterns, step counts, sleep duration, phone usage statistics, wearable HRV, PHQ-8 baseline, and (when available) voice prosody and keyboard linguistic features.

A higher score indicates behavioral patterns that have been statistically associated with elevated mental health risk in population-level research data.

This is NOT a diagnosis. The score is an awareness indicator. Significant or sustained score changes warrant discussion with a mental health professional.

6.3 Risk Trend Chart

The 7-day risk trend line chart shows how your composite risk score has changed over the past week. Each data point represents the risk score calculated for that day.

- Days with data: filled circle markers on the teal line
- Today: highlighted with a label
- A status pill (e.g., 'Risk down vs yesterday') appears in the top-right corner

6.4 Today's Highlights

Four small insight cards provide at-a-glance comparisons:

Card	What It Shows
Improvement compared to yesterday	Whether today's risk is lower or higher than yesterday's
Better / Worse than weekly average	Comparison of today's risk to your 7-day average
Mood lifting recently	3-day mood average vs current mood (e.g. '3-day avg 3.1 vs 2.8')
Activity status	Step goal progress for today (e.g. 'Goal 10000 steps · today 0')

6.5 Scuba Recommendations

At the bottom of the Home screen, ScubaMind displays 2–3 personalized, actionable recommendations derived from your current behavioral data. Examples include:

- 'Gentle movement can ease low mood — Your activity is low, and risk signals are elevated. Even a 10-minute walk outside activates mood-lifting neurotransmitters.'

- 'Write a 2-minute journal — Note one thing that went well today and one thing you're grateful for. Even a brief reflection can meaningfully shift your mood.'

Recommendations update automatically as new data is collected and analyzed.

6.6 Logging Your Mood

Tap Log mood on the Home screen to manually record how you are feeling right now. You will be taken to the Mood logging screen, where you can select from five emotional states represented as ScubaMind mascot emoji faces.

- Face 1 (leftmost): Very low / Struggling
- Face 2: Low / Difficult
- Face 3 (center): Neutral / Okay
- Face 4: Good / Positive
- Face 5 (rightmost): Great / Excellent

Your mood entries are stored locally and contribute to the mood signal on the Home screen and Mood detail page.

6.7 Pull-to-Refresh

To manually refresh all data on the Home screen, pull down on the screen. This triggers a fresh load of all behavioral data modules, including steps, sleep, and risk score. A brief loading indicator will appear on the stat cards as data refreshes.

7 Goals

The Goals screen helps users create, follow, and maintain personalized wellness objectives as part of their daily mental well-being routine. From this page, users can view their active goals, check current progress, update completed actions, create new goals, explore suggested goals, manage voice check-ins, and review previously completed goals.

ScubaMind supports different goal categories such as activity, sleep, mindfulness, nutrition, hydration, and voice journaling. Users can also enable daily or weekly reminders so that the application sends notifications at the selected time. When a goal is completed, ScubaMind may reward the user with an underwater-themed achievement badge, helping users stay motivated and consistent over time. Tap the Goals icon in the bottom navigation bar to access this section.

The Goals screen is designed not only for task tracking but also for encouraging long-term consistency and positive behavior change. Instead of presenting goals as simple checklist items, ScubaMind connects them with progress indicators, reminders, history tracking, and achievement badges. This helps users understand how small actions, such as drinking water, improving sleep routine, completing a short voice check-in, or taking a walk, can contribute to their overall well-being patterns. Completed goals are stored in the user's history, and recurring goals can support daily or weekly routines without requiring users to recreate the same objective repeatedly.

ScubaMind



My Well-Being Goals

Active Goals Edit

Daily Voice Check
0 of 30 sec recorded

0% >

Drink Water
0 of 8 glasses done

0% >

Eat Vegetables
0 of 5 servings done

0% >

Sleep Routine
0 of 8 hours done

0% >

My Voice Recordings

>

PHQ-8 Check-in

>



Eat Vegetables
0 of 5 servings done

0% >

Sleep Routine
0 of 8 hours done

0% >

My Voice Recordings

>

PHQ-8 Check-in
Quick mood questionnaire

>

Discover New Goals

Physical

Mental

Rest

Nutr

Sleep Routine
8 hours every night

Add to Goals

Screen-Free Time
1 hour before bed

Add to Goals



← Goal History

Your Journey
Track completed goals and voice progress in one place.

1 Completed

1 Active

0 Missed

2 Total

2 entries

This overview

Completed

Active

Missed

1 • 50%

1 • 50%

0 • 0%

Completed Goals

Voice Recordings

PHQ Results

Reactivate Goals

Personal

Daily Walking
30 min completed

Reactivate





NEW BADGE EARNED!



HYDRATION GOAL
COMPLETE #1

Great job! Your hydration goal is completed.

Tap to continue

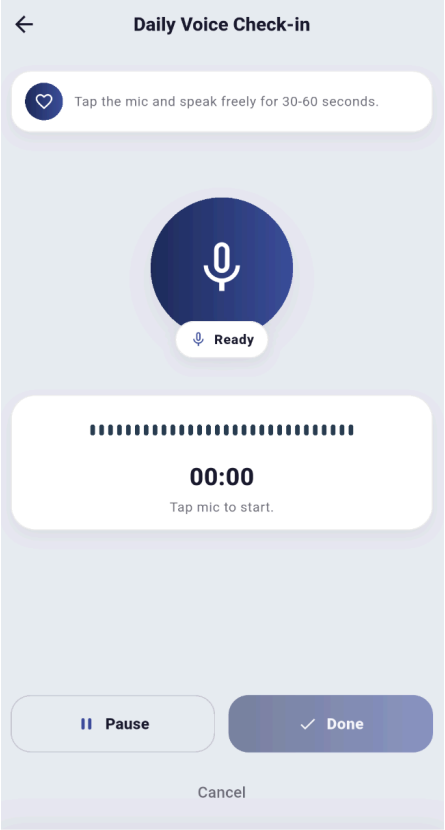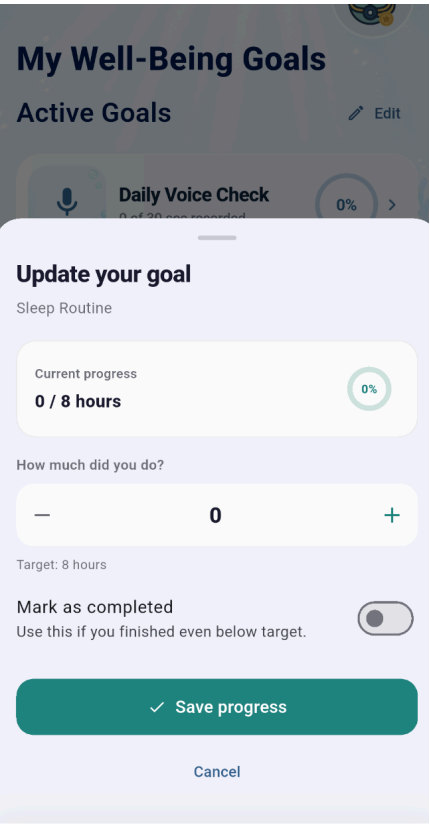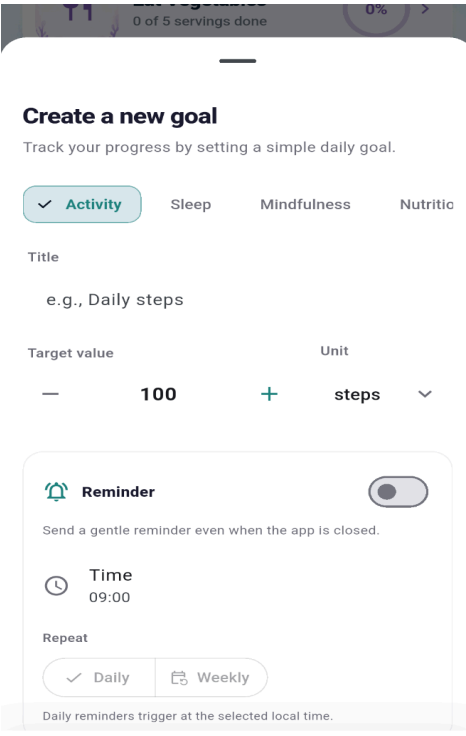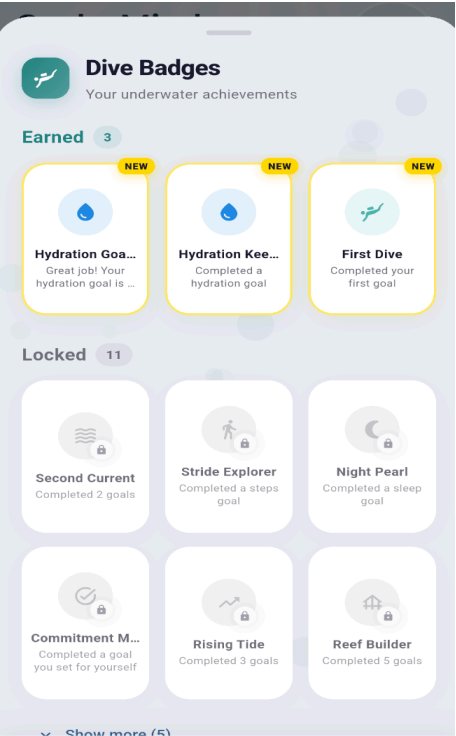

7.1 Active Goals

The top section of the Goals screen shows your currently active goals. Each goal card displays:

- Goal icon and title (e.g., 'Daily Voice Check')
- Current progress (e.g., '0 of 30 sec recorded')
- A circular progress indicator showing completion percentage
- A '>' arrow to open the goal detail/update screen

Default pre-configured goals include:

Goal	Default Target
Daily Voice Check	30 seconds of voice recorded
Drink Water	8 glasses per day
Eat Vegetables	5 servings per day
Sleep Routine	8 hours of sleep

7.2 Updating Goal Progress

To record progress on a goal:

21. Tap the '>' arrow on any active goal card, or tap the goal title.
22. The 'Update your goal' sheet slides up from the bottom.
23. You will see your Current progress (e.g., '0 / 8 hours') and a percentage ring.
24. Use the – and + buttons to enter how much you completed (e.g., 6 hours of sleep).
25. To mark as fully completed regardless of the number, toggle Mark as completed.
26. Tap Save progress. The goal card will update immediately.
27. Tap Cancel to discard changes.

7.3 Creating a New Goal

Tap the green + button in the bottom-right corner of the Goals screen to open the 'Create a new goal' sheet.

Field	Description
Category	Choose from Activity, Sleep, Mindfulness or Nutrition (tab selector at the top)
Title	Enter a name for your goal (e.g., 'Daily steps')
Target value	Use – and + buttons to set the numeric target (e.g., 100)
Unit	Select the unit from the dropdown (steps, hours, glasses, servings, minutes, etc.)
Reminder (toggle)	Enable receiving a notification when it is time to check in on this goal

Field	Description
Reminder Time	Set the time of day for the reminder (default: 09:00)
Repeat	Choose Daily or Weekly repetition for the reminder

After filling in all fields, scroll down and tap Save Goal (or equivalent confirm button) to add it to your active goals.

7.4 Discover New Goals

Scroll down on the Goals screen to reach the 'Discover New Goals' section. Goals are organized by category tabs:

- Physical — exercise, steps, movement-related goals
- Mental — mindfulness, journaling, relaxation goals
- Rest — sleep and screen-free time goals
- Nutrition — hydration, diet-related goals

Tap Add to Goals on any suggestion card to instantly add that goal to your active goals list.

7.5 Voice Recordings

Access from Goals → My Voice Recordings. This screen shows all your saved voice journal entries.

Each recording card shows:

- Recording title and date
- Duration and file size
- Transcription status: a loading spinner with 'Transcribing' while processing, or the transcribed text (up to 5 lines) when complete

Actions available on the recording list:

- Tap Play to listen to the recording — icon toggles to Pause while playing
- Tap Pause to pause playback
- Tap the Delete (bin) icon → confirm in the dialog to permanently delete a recording.

Transcription Processing Time

The on-device Whisper ASR model processes each 30-second clip in approximately 5–10 seconds on a mid-range device.

During processing, the recording card shows a 'Transcribing' indicator. Do not delete the recording during this time.

If a recording is silent, the transcription will return an empty string (this is normal).

7.6 PHQ-8 Check-in

A 'PHQ-8 Check-in — Quick mood questionnaire' entry is available from the Goals screen. Tap it to complete a new PHQ questionnaire at any time. Results are stored locally and contribute to the AI model's calibration.

7.7 Daily Voice Check-in

Tap the Daily Voice Check goal or navigate to it to reach the Voice Check-in recording screen:

- Tap the large microphone button (shows 'Ready') to begin recording.
- Speak freely for 30–60 seconds. The timer and audio waveform are visible.
- Tap Pause at any time to pause; tap the button again to resume.
- Tap Done when finished to save the recording.
- Tap Cancel to discard without saving.

After saving, ScubaMind processes the audio on-device using the Whisper ASR model and generates a text transcription. This may take up to 10 seconds. The transcription appears on the recording card once complete.

7.8 Goal History & Badges

Tap the history/edit icon in the Goals screen to view Goal History. This screen shows:

Section	Description
Summary Tabs	Filter by Completed Goals, Voice Recordings, or PHQ Results
Statistics Ring	Donut chart showing Completed / Active / Missed proportions
Reactivate Goals	List of previously completed goals with a Reactivate button
Dive Badges	Your earned achievements — see below

7.9 Dive Badges — Achievements

ScubaMind uses an underwater-themed achievement system to recognize consistent healthy behaviors:

Badge	How to Earn
First Dive	Complete your very first goal
Hydration Goal Complete #1	Complete a hydration (water drinking) goal
Hydration Keeper	Complete multiple hydration goals
Stride Explorer	Complete a step / activity goal
Night Pearl	Complete a sleep goal
Second Current	Complete 2 goals total
Rising Tide	Complete 3 goals total
Reef Builder	Complete 5 goals total
Commitment Master	Complete a goal you set for yourself

When you earn a new badge, a full-screen celebration animation appears ('NEW BADGE EARNED!') showing the badge icon and its name. Tap anywhere to continue.


All badges are visible in the Dive Badges sheet (accessible from Goal History). Earned badges are highlighted with a 'NEW' tag; locked badges are shown greyed out with a padlock.

8 Chat - Scuby AI Companion

The Chat screen provides access to Scuby, ScubaMind's AI-powered conversational companion. Scuby is powered by the Google Gemini API and is designed to offer empathetic, reflective, and supportive dialogue. Through this screen, users can type messages, ask for wellbeing suggestions, use quick action buttons such as sleep tips or breathing support, and receive personalized responses based on high-level behavioral summaries. Scuby is intended to help users reflect on their current mood, understand possible lifestyle patterns, and take small supportive actions such as setting a goal, journaling briefly, or taking a short break.

8.1 Chat Screen Layout

Element	Description
Avatar Area	Avatar Area — A 3D animated character (“Scuby”) — an AI companion avatar with warm and supportive aesthetics.
Online Status	Green dot with 'online' label
Top-right Icons	Video (camera), Chart (behavioral data), Delete (clear conversation) icons
Quick Action Pills	Sleep tips / Set a goal / Breathing — one-tap conversation starters
Message List	Chronological conversation history (Scuby on left, your messages on right in teal)
Input Bar	+ icon / Text field ('Write to SCUBA Assistant...') / Microphone / Send button

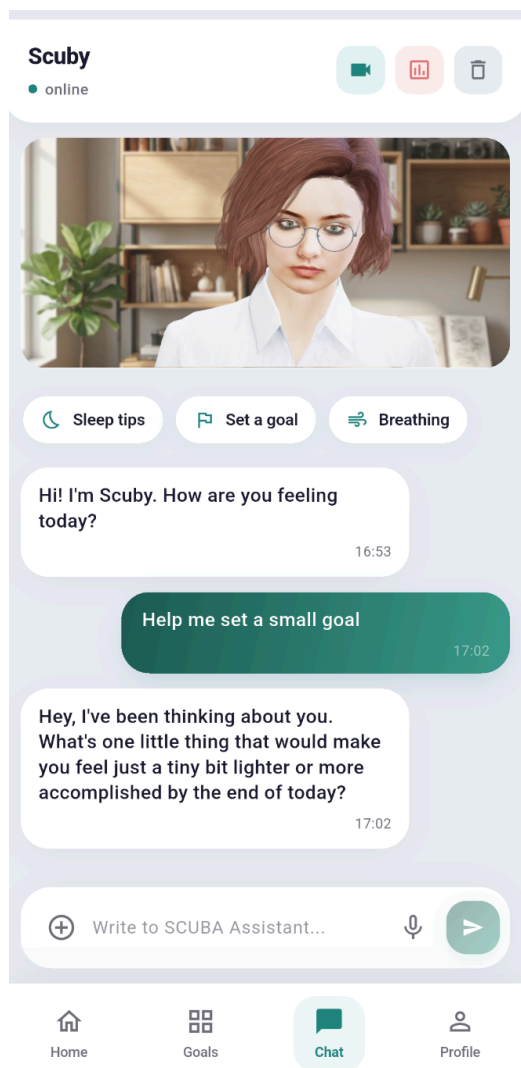

8.2 Starting a Conversation

Scuby greets you automatically: 'Hi! I'm Scuby. How are you feeling today?'

You can interact in three ways:

28. Type a message in the text field and tap the teal Send (arrow) button.
29. Tap the microphone icon to send a voice message.
30. Tap one of the Quick Action pills (Sleep tips, Set a goal, Breathing) to start a predefined conversation flow.

8.3 How Scuby Uses Your Data

Before each conversation turn, ScubaMind automatically constructs a personalized context prompt by aggregating recent behavioral insights from your local database. This context includes high-level summaries such as:

- Your recent activity levels compared to your weekly average
- Your sleep duration over the past 3 days
- Your current risk level and any notable trend changes

This allows Scuby to give you relevant, personalized responses — for example, if your activity has been low, Scuby may suggest a short walk or ask about your energy levels.

What Is and Is Not Sent to Gemini API

SENT: High-level, anonymized behavioral summaries (e.g., 'activity level: low this week').

NEVER SENT: Raw voice recordings, transcription text, keyboard text, exact GPS coordinates, or any personally identifiable information.

Only the messages you type directly into the Scuby Chat and high-level anonymized behavioral summaries may be sent to the Gemini API to generate a response. Raw voice recordings, exact GPS coordinates, keyboard text from other apps, full sensor logs, and identifiable health records are not sent.

8.4 Safety Guardrails

All responses from Scuby are processed through a SafetyGuardrailLayer that:

- Monitors for high-risk language indicating self-harm or crisis
- Immediately redirects to emergency resources if such content is detected
- Uses Gemini's built-in safety filters in combination with custom keyword monitoring

⚠ If You Are in Crisis

Scuby is not a crisis service. If you are in immediate danger or experiencing thoughts of self-harm, contact local emergency services immediately.

Turkey: 112 Emergency Call Center | International: findahelpline.com






Scuby will provide crisis resource information if it detects high-risk content.

8.5 Clearing a Conversation

Tap the Delete (bin) icon in the top-right corner of the Chat screen to clear your current conversation history. This action is irreversible.

9 Profile

The Profile screen brings together the user's personal progress, wellbeing status, and application controls in one place. It allows users to review their level, XP progress, streak count, check-ins, earned badges, and active goals while also providing quick access to mood logging, recommendations, privacy settings, app settings, and support options. This page helps users understand their overall engagement with ScubaMind and makes it easier to manage both personal well-being features and technical preferences without navigating through multiple screens.

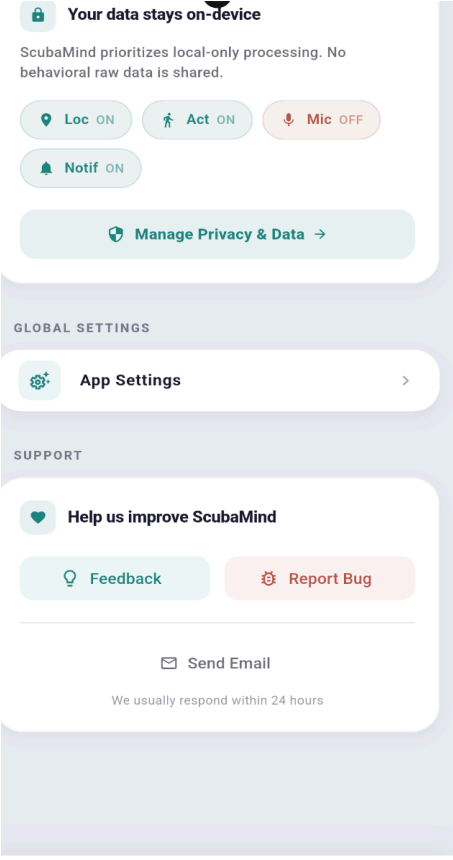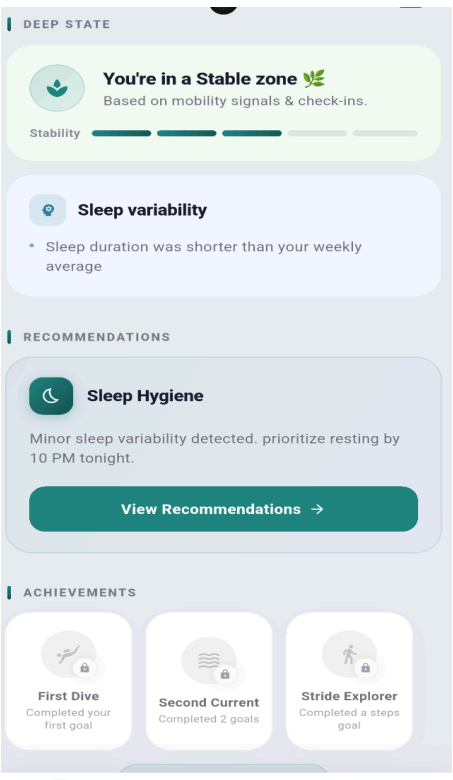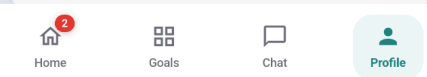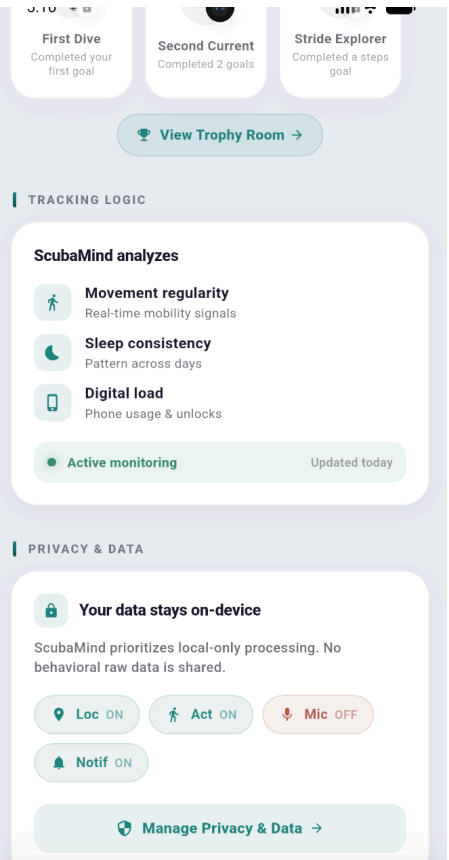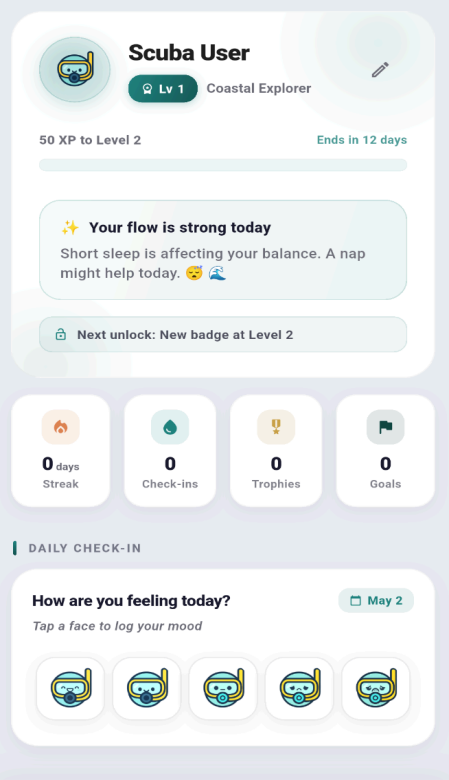

9.1 User Card

Element	Description
Avatar	Your ScubaMind mascot avatar (scuba diver character). Tap the pencil icon to edit your profile.
Username	Default: 'Scuba User'
Level Badge	Current level (e.g., 'Lv 1 Coastal Explorer') with XP progress bar
XP Progress	'50 XP to Level 2' with a progress bar and 'Ends in N days' countdown
Daily Insight	A short personalized message (e.g., 'Your flow is strong today — Short sleep is affecting your balance. A nap might help today.')
Next Unlock	Shows what reward you will unlock at the next level

9.2 Stats Bar

Metric	Description
Streak (days)	Number of consecutive days you have completed at least one goal
Check-ins	Total number of PHQ / mood check-ins completed
Trophies	Total number of badges earned
Goals	Number of active goals currently set

9.3 Daily Check-in

The Daily Check-in section appears below the stats bar. It shows today's date and five mascot faces representing mood levels.

Tap a face to log your mood for today:

- Face 1 (Sad/Struggling): Very low mood
- Face 2 (Slightly sad): Below average
- Face 3 (Neutral): Neither good nor bad
- Face 4 (Happy): Positive mood
- Face 5 (Very happy): Excellent mood

9.4 Deep State

The Deep State section provides an overview of your current behavioral stability:

- Stability Zone indicator (e.g., 'You're in a Stable zone') based on mobility signals and check-ins
- Stability bar showing your current position on a stability spectrum
- Sleep Variability indicator — flags if your recent sleep duration deviated from your weekly average

9.5 Recommendations

A personalized recommendation card appears in the Recommendations section. For example: 'Sleep Hygiene — Minor sleep variability detected. Prioritize resting by 10 PM tonight.'

Tap View Recommendations → to see all current recommendations.

9.6 Achievements

A preview of your earned and locked badges is shown. Tap View Trophy Room → to open the full Dive Badges sheet with all earned and locked achievements (see Chapter 7.9).

9.7 Tracking Logic

The Tracking Logic section shows which data sources are currently active:

- Movement regularity — Real-time mobility signals
- Sleep consistency — Pattern across days
- Digital load — Phone usage and unlocks
- Active monitoring — Updated today (or last update date)

9.8 Privacy & Data (Quick Access)

The Privacy & Data section in the Profile screen shows your current permission status at a glance:

- Loc ON/OFF — Location tracking
- Act ON/OFF — Activity (step) tracking
- Mic ON/OFF — Microphone/voice recording
- Notif ON/OFF — Notifications

Tap Manage Privacy & Data → to open the full Privacy & Data settings screen (see Chapter 12).

9.9 Global Settings

Tap App Settings → to open the Settings screen (see Chapter 12).




9.10 Support

The Support section at the bottom of the Profile screen provides:

- Feedback — Submit general feedback about the application
- Report Bug — Report a technical issue with the app
- Send Email — Contact the development team directly (response within 24 hours)

10 Risk Breakdown

The Risk Breakdown screen provides a detailed analysis of the factors contributing to your overall behavioral risk score. Access it by tapping 'View breakdown >' on the Home screen's Mental Risk Level card.

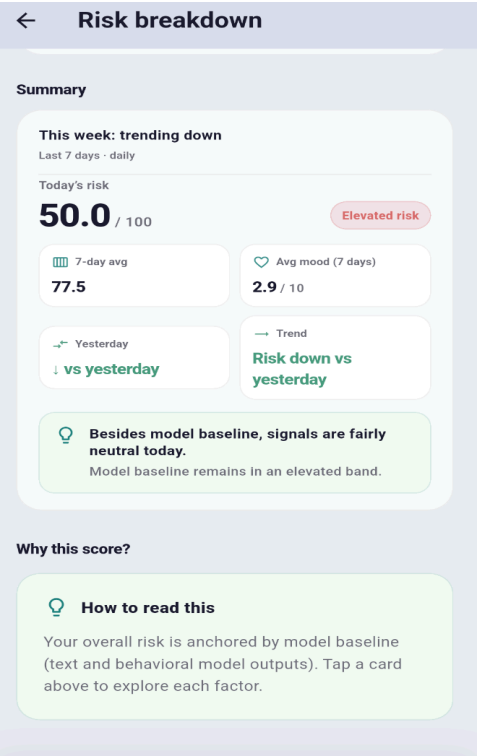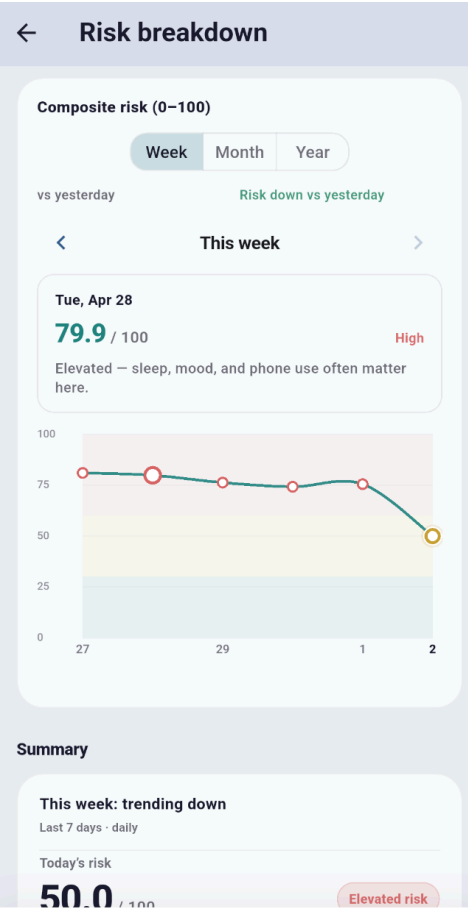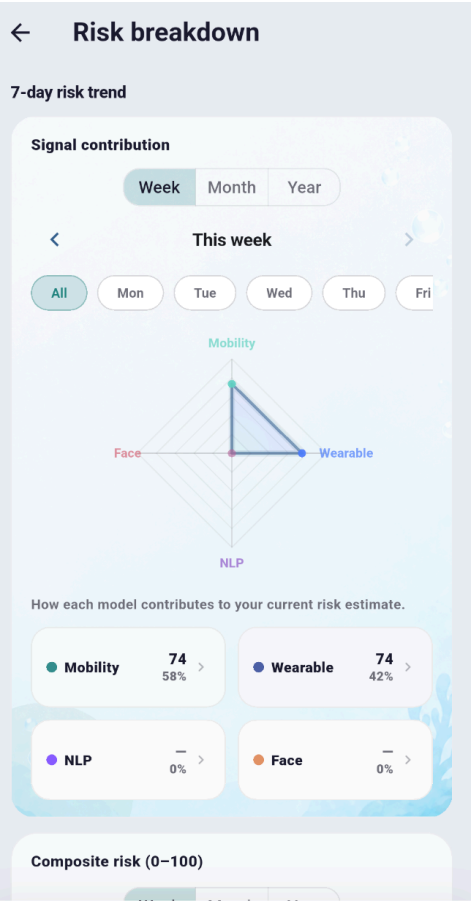

10.1 Signal Contribution Radar

The top section shows a radar (spider) chart with four axes:

- Mobility — Contribution from GPS-based movement patterns
- Wearable — Contribution from steps and sleep data
- NLP — Contribution from text analysis (keyboard/voice transcription)
- Face — Contribution from facial expression analysis (if enabled)

Use the Week / Month / Year toggle to change the time window. Use < This week > navigation arrows to view previous periods. You can also filter by individual day (Mon, Tue, Wed, Thu, Fri tabs) or view all days combined.

Below the radar chart, four signal cards show:

- Mobility: score/contribution percentage
- Wearable: score/contribution percentage
- NLP: score/contribution percentage
- Face: score/contribution percentage

Tap any signal card to navigate to the detailed signal page for that modality.

10.2 Composite Risk Trend

The lower section of the Risk Breakdown screen shows the Composite Risk (0–100) trend chart with the same Week / Month / Year view options.

Tapping on a specific day's data point shows a detail card with:

- Date (e.g. Tue, Apr 28)
- Risk score (e.g. 79.9 / 100)
- Qualitative label (e.g. 'High')
- Brief explanation (e.g., 'Elevated — sleep, mood, and phone use often matter here.')

10.3 Summary Section

Field	Description
Weekly Trend	Direction of change (e.g., 'This week: trending down')
7-day Average	Mean risk score over the past 7 days
Average Mood (7 days)	Mean mood rating over the past 7 days (scale 1–10)
vs Yesterday	Whether today's risk is higher or lower than yesterday's
Trend	Direction badge (e.g., 'Risk down vs yesterday')
Status Summary	Qualitative note (e.g., 'Besides model baseline, signals are fairly neutral today.')

10.4 Why This Score? — How to Read This

At the bottom of the Risk Breakdown page, the 'Why this score?' section explains the logic behind your risk score:

'Your overall risk is anchored by the model baseline (text and behavioral model outputs). Tap a card above to explore each factor.'


Each modality card links to its dedicated signal page, where you can examine the underlying raw data and trends.

11 Detailed Signal Pages

Each data modality collected by ScubaMind has its own detailed view. Access these pages by tapping the relevant signal card on the Risk Breakdown screen or through the mini-stat cards on the Home screen.

11.1 Activity Page

Shows your step count and movement data.


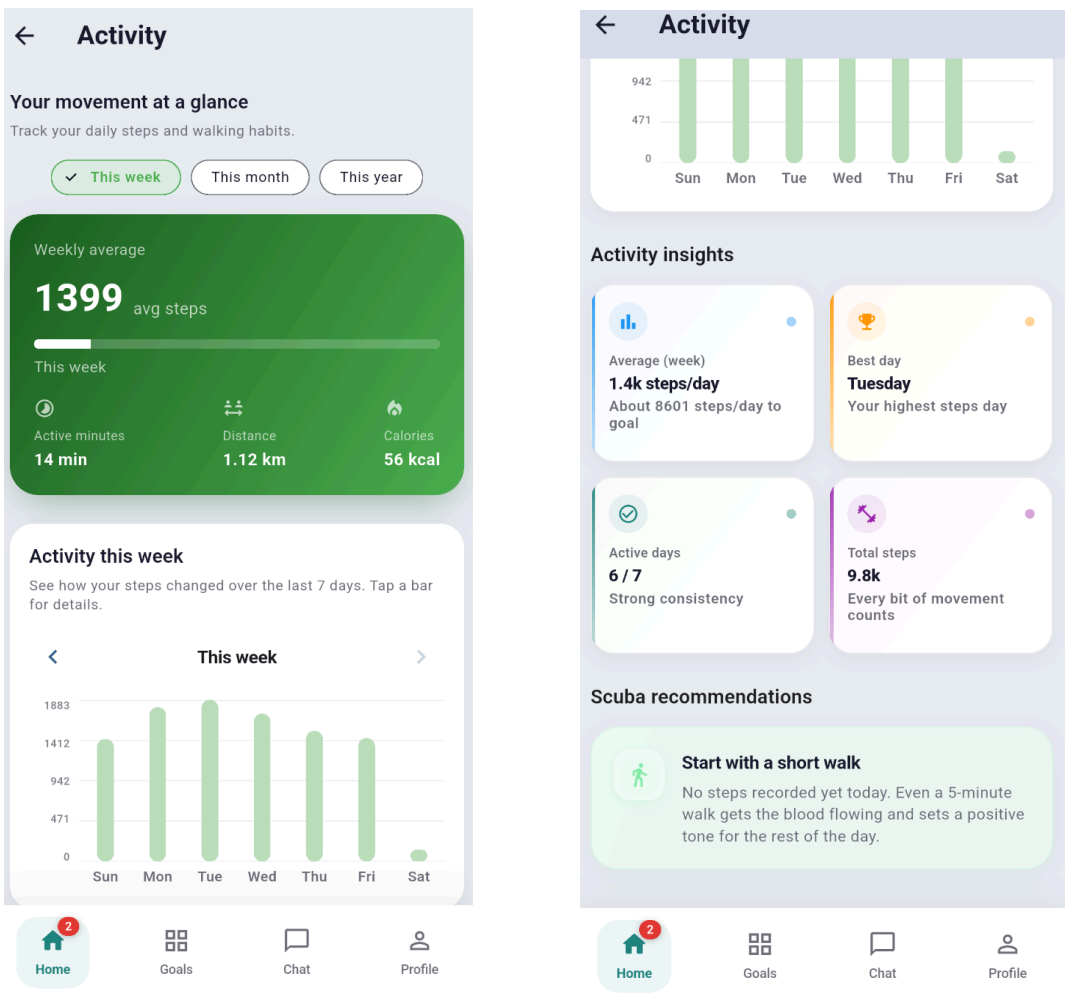


Section	What It Shows
Weekly Average Card	Average daily steps, active minutes, distance, and calories for the selected period
Activity This Week Chart	Bar chart of daily step counts. Tap a bar to see that day's detail.
Activity Insights	Four cards: average steps/day, best day of the week, active days (X/7), total steps
Scuba Recommendations	Contextual activity recommendations (e.g., 'Start with a short walk')

Use the This week / This month / This year toggles and < > navigation arrows to explore different time periods.

11.2 Sleep Page

Shows your sleep duration and patterns.

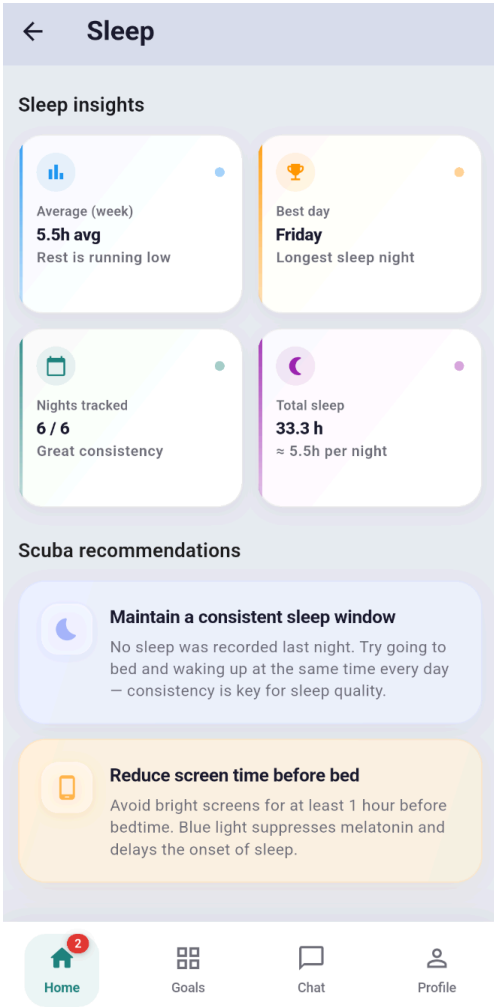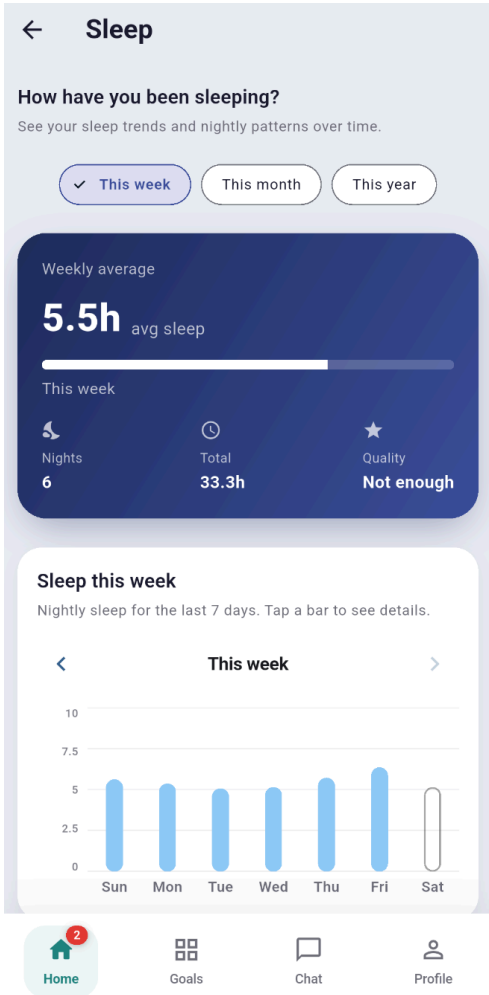


Section	What It Shows
Weekly Summary Card	Average sleep duration, number of nights tracked, total sleep hours, quality rating

Section	What It Shows
Sleep This Week Chart	Bar chart of nightly sleep duration. Tap a bar for that night's detail.
Sleep Insights	Average weekly sleep, best sleep night (day with most hours), nights tracked, total sleep
Scuba Recommendations	Sleep hygiene tips (e.g., 'Reduce screen time before bed')

Sleep data is inferred from screen-on/off events and, if connected, from Samsung Galaxy Watch sleep stage data via Health Connect.

11.3 Mood Page

Shows your mood log history and emotional patterns.

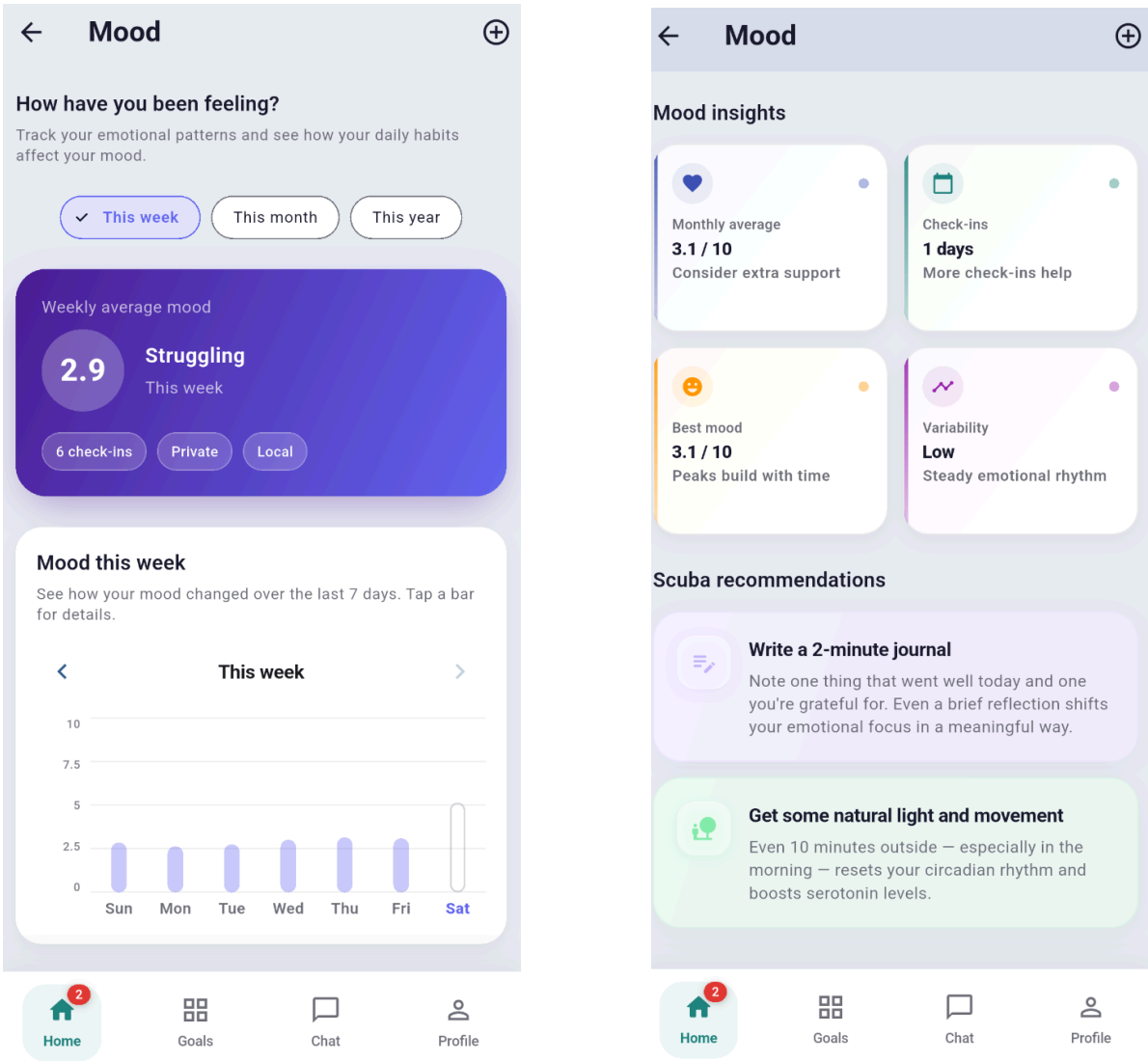


Section	What It Shows
Weekly Average Card	Mean mood score for the period (scale 1–10), qualitative label (e.g. 'Struggling'), number of check-ins
Mood This Week Chart	Bar chart of daily mood values. Tap a bar for that day's mood.

Section	What It Shows
Mood Insights	Four cards: monthly average, check-in count, best mood value, variability level
Scuba Recommendations	Mood-specific suggestions (e.g., 'Write a 2-minute journal')

11.4 Phone Usage Page


Shows your screen time and digital behavior patterns.

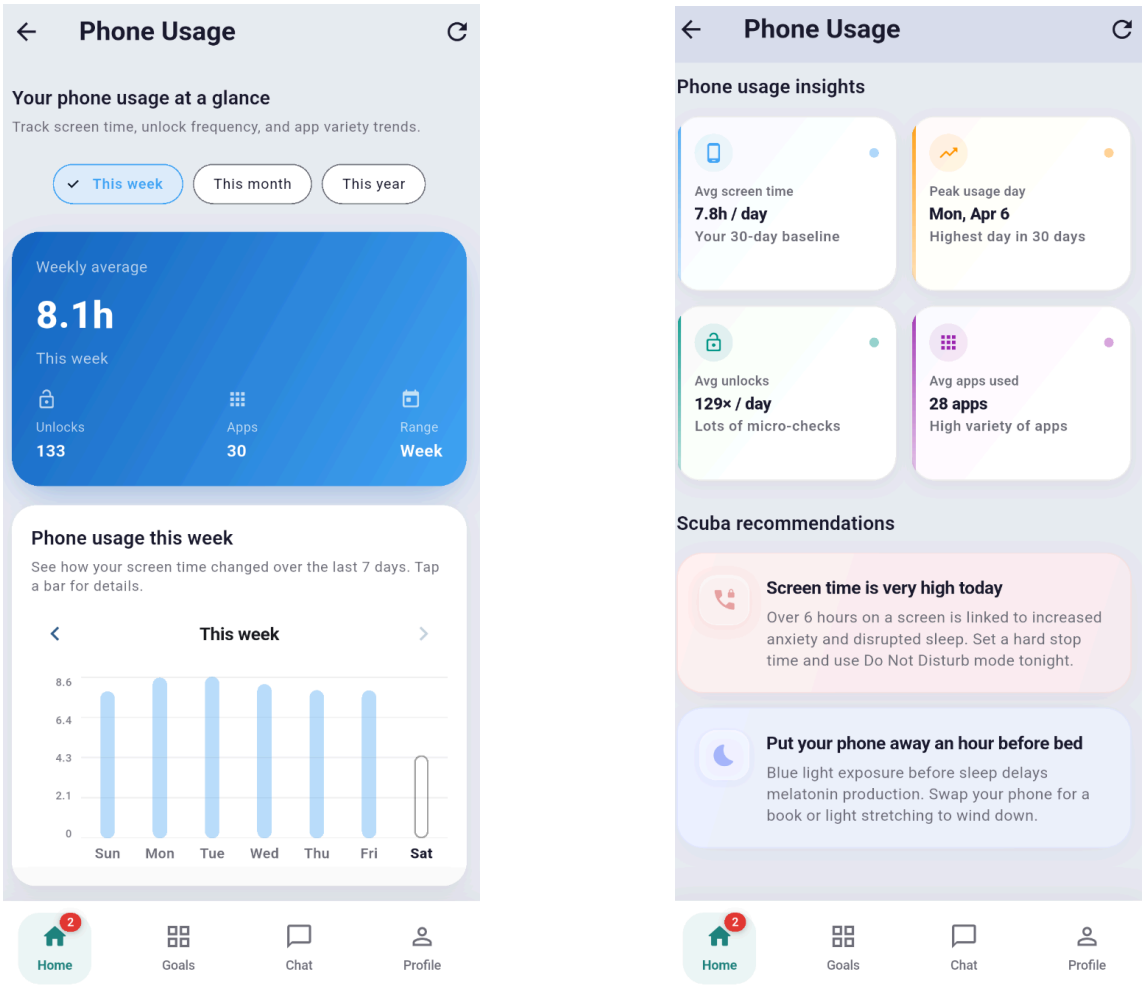

Section	What It Shows
Weekly Average Card	Average daily screen time in hours, total unlocks, apps used, time range
Phone Usage This Week Chart	Bar chart of daily screen time. Tap a bar for that day's detail.
Phone Usage Insights	Average screen time, peak usage day, average daily unlocks, average apps used
Scuba Recommendations	Digital wellbeing tips (e.g., 'Screen time is very high today — Set a hard stop time')

11.5 Mobility Page

Shows your GPS-based movement patterns and location diversity.

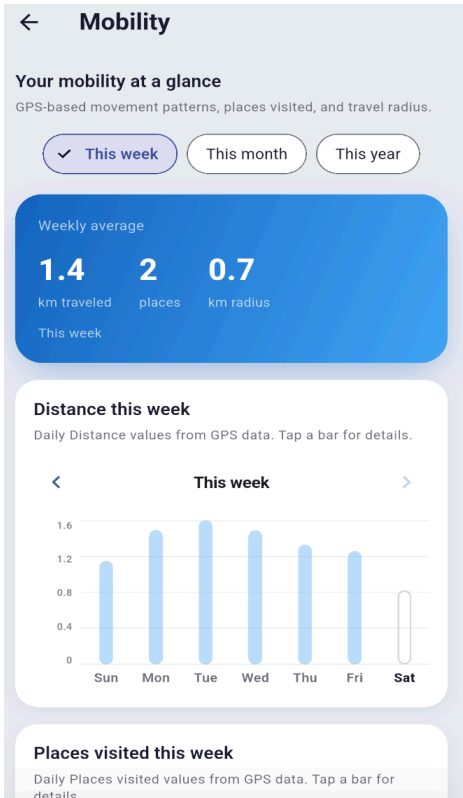

Section	What It Shows
Weekly Summary Card	Average km traveled, number of distinct places visited, mobility radius (km)
Distance This Week Chart	Bar chart of daily distance in km. Tap a bar for that day.
Places Visited Chart	Bar chart of distinct places visited per day.
Mobility Radius Chart	Bar chart of daily radius of gyration values.
Top Places (7 days)	Privacy-friendly location clusters (~100m radius). Most frequent location labeled ' Likely home'.
Mobility Insights	Daily average distance, most active day, place diversity insights

Mobility Privacy Protection

Top places are shown as approximate coordinates clustered to ~100m radius. Exact addresses are never stored.


Location data never leaves your device.

You can disable Location tracking at any time from Privacy & Data.

11.6 Mobility Signal (AI Analysis) Page

This page, accessible from the Risk Breakdown screen, shows the AI model's interpretation of your mobility data:



- Fallback estimate tag (shown when there is insufficient data for a full model estimate)
- Mobility risk estimate (e.g., 82 / 100)
- Routine average, Places average, Distance average, PHQ estimate
- Mobility Risk Trend chart
- Routine Stability chart over time
- Mobility Composition breakdown (Routine stability, Place diversity, Radius of gyration, Home range, Movement variability)
- Mobility Intensity chart color-coded Low / Medium / High per day)
- Key Indicators: Routine stability, Place diversity, Radius of gyration, Home range, Total distance, Significant places

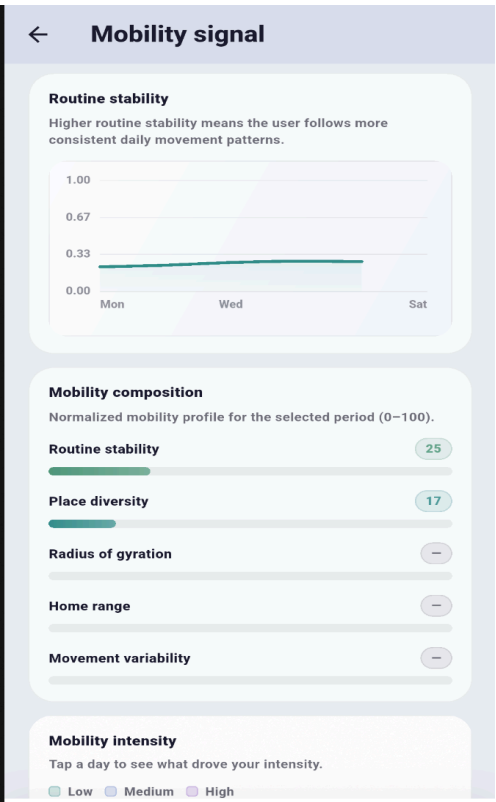

11.7 Wearable Signal Page


Shows data imported from Samsung Galaxy Watch via Health Connect.

Section	What It Shows
Wearable Estimate Card	AI risk estimate from wearable data, contribution %, steps avg, sleep avg, data coverage
Activity + Recovery Chart	Dual-axis bar (steps) and line (sleep hours) overlay chart
Daily Rings	Three concentric ring chart: Activity (steps), Recovery (sleep duration), Consistency (day-to-day stability)
Weekly Activity Strip	Day-by-day tabs (Mon–Sat). Tap a day to see that day's steps and sleep.
Key Indicators	Activity regularity, Recovery support level, Consistency, Confidence

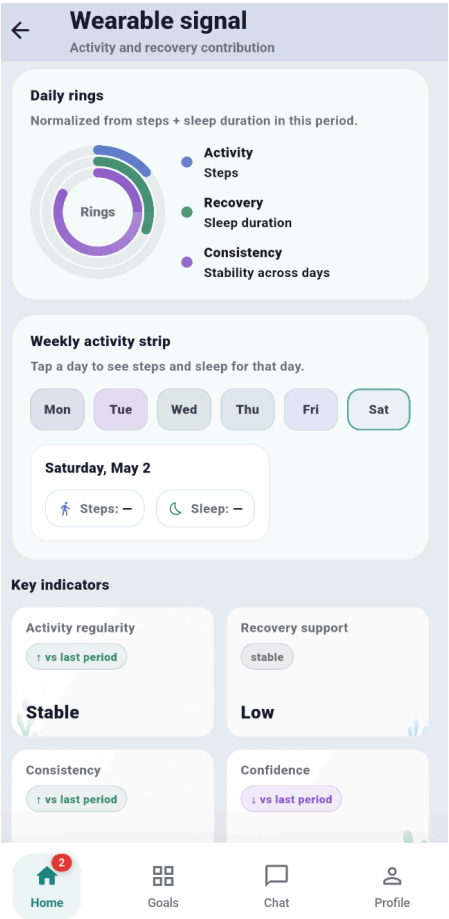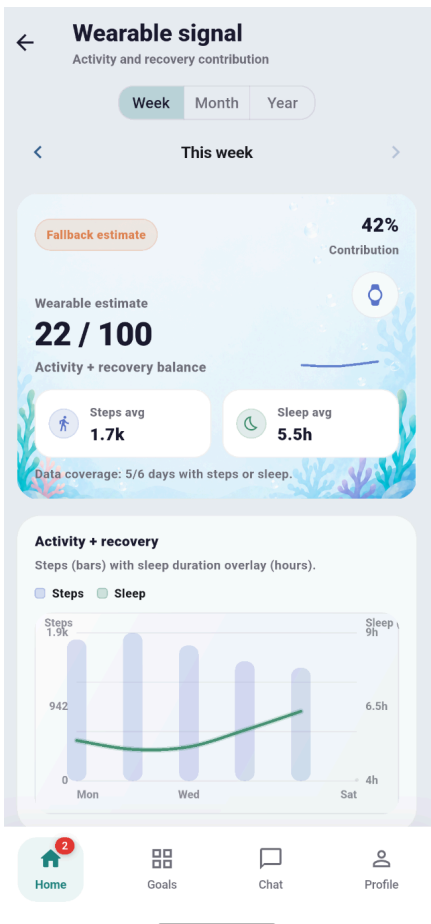

12 Privacy & Data

The Privacy & Data screen gives you full control over what data ScubaMind collects and how it is used. Access it from Profile → Privacy & Data → Manage Privacy & Data.

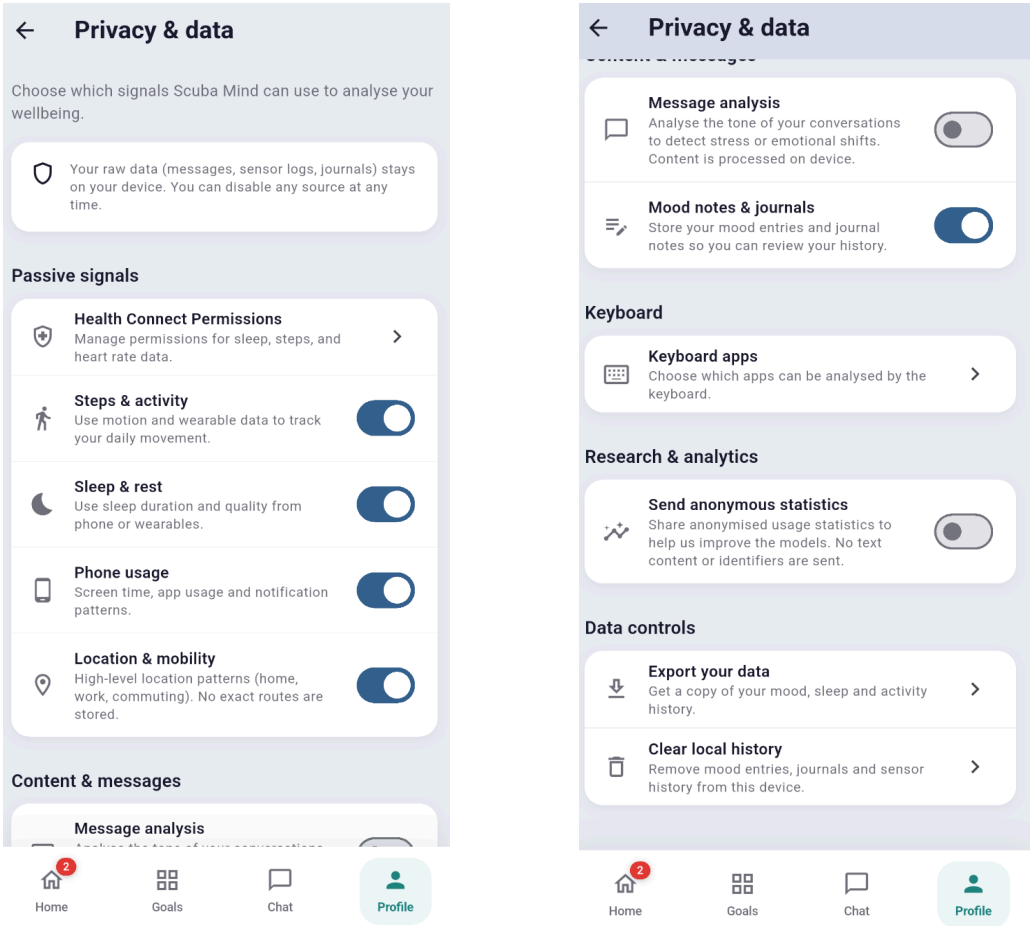

12.1 Passive Signals

Toggle each passive data source on or off individually:

Signal	Description
Health Connect Permissions	Tap to manage sleep, steps, and heart rate permissions in the Health Connect system settings
Steps & activity	Use motion and wearable data to track your daily movement
Sleep & rest	Use sleep duration and quality estimates from phone or wearables
Phone usage	Monitor screen time, app usage patterns, and notification frequency
Location & mobility	High-level location patterns (home, work, commuting). No exact routes are stored.

12.2 Content & Messages

Feature	Description
Message analysis	Analyze the tone of your text conversations to detect stress or emotional shifts. All content processed on-device.
Mood notes & journals	Store your mood entries and journal notes so you can review your history.

12.3 Keyboard

Tap Keyboard apps → to choose which apps are analyzed by the ScubaMind custom keyboard. You can add or remove specific apps from the monitored list at any time.

12.4 Research & Analytics

Feature	Description
Send anonymous statistics	Share anonymized, aggregated usage statistics to help improve the AI models. No text content or identifiers are included in this sharing. Toggle OFF by default.

12.5 Data Controls

Option	What It Does
Export your data	Generate a PDF wellbeing report for a chosen time period (see Chapter 15)
Clear local history	Remove mood entries, journals, and sensor history from your device

12.6 Complete Data Purge (Delete All Data)

Navigate to Profile → Privacy & Data → Clear local history → Delete All Data for a complete and irreversible data purge.

This action is irreversible

Deleting all data performs a cryptographic key destruction operation on the AES-256 encrypted database.





All stored behavioral data, PHQ responses, voice recording metadata, goals, and mood entries are permanently destroyed.

The encryption keys are removed from the Android Keystore, making data unrecoverable.

The app will restart in a fresh state, equivalent to a first-time installation.

13 Settings

Access Settings from Profile → App Settings.

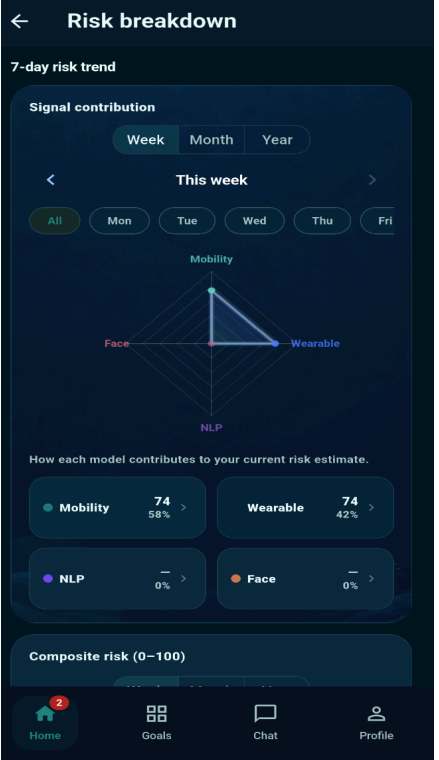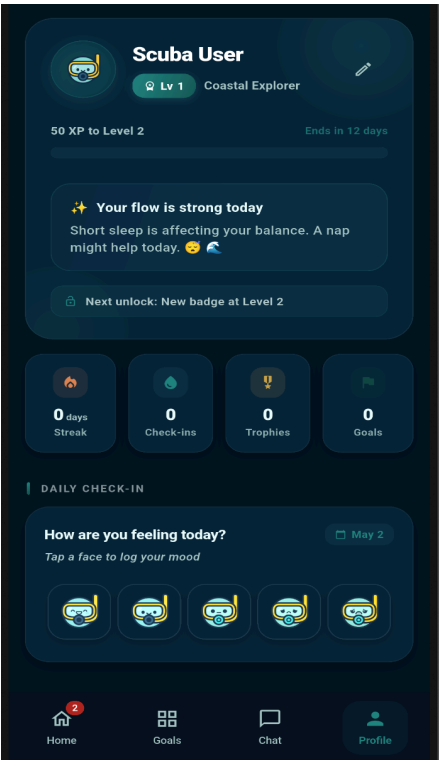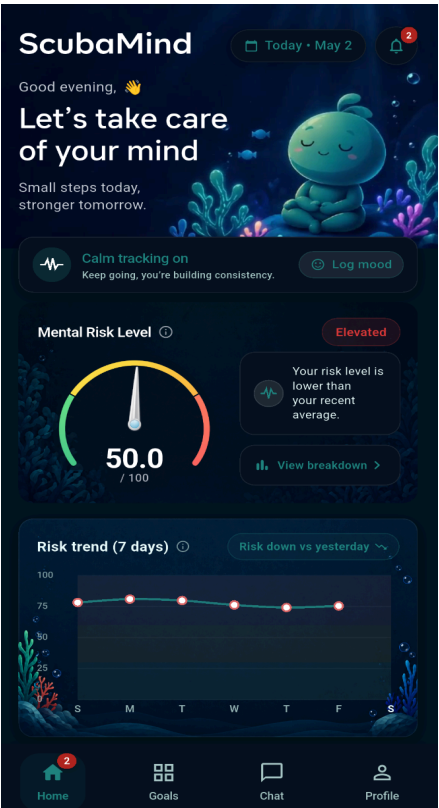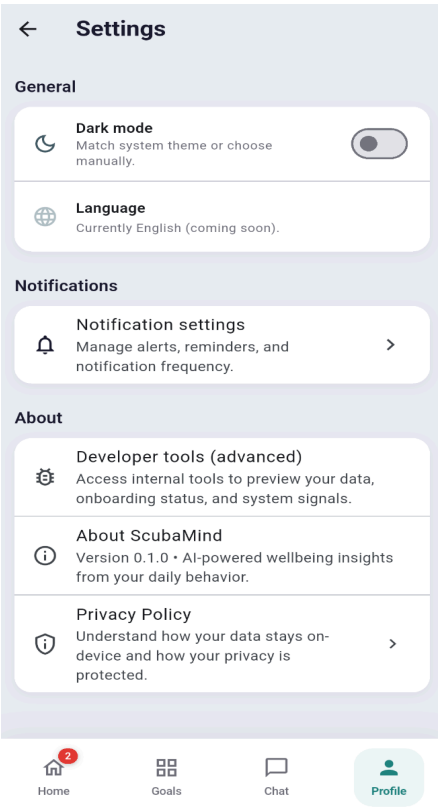

13.1 General

These settings allow users to personalize their experience without changing the core functionality of the application. For example, users can switch between light and dark appearance modes depending on their visual preference or device theme, and they can review the currently supported language options. This section is intended to make the app more comfortable, accessible, and user-friendly for everyday use.

The Dark mode option is especially useful for users who prefer a lower-brightness interface, particularly in the evening or in low-light environments. When enabled, ScubaMind adjusts its colors and cards to a darker visual theme while preserving readability and contrast. Users may either let the app follow the system theme automatically or manually choose the display mode they prefer. This helps provide a more flexible and visually comfortable experience across different usage conditions.

Setting	Description
Dark mode	Toggle between light and dark themes. Options: Match system theme, or enable manually. All screens in ScubaMind support dark mode.
Language	Currently, English only. Additional language support (including Turkish) is planned for a future release.

13.2 Notifications

Tap Notification settings → to manage:


- Goal reminder alerts and their frequency
- Daily check-in prompts
- System status notifications from background services

13.3 About

Item	Description
Developer tools (advanced)	Access internal diagnostic tools to preview your raw data, onboarding status, and system signals. Intended for development and debugging only.
About ScubaMind	App version (0.1.0), build info, and feature tagline
Privacy Policy	Read the full privacy policy explaining how your data stays on-device and how your privacy is protected.

14 Notifications

The Notifications page displays all in-app alerts and reminders. Access it by tapping the bell icon in the top-right corner of the Home screen.

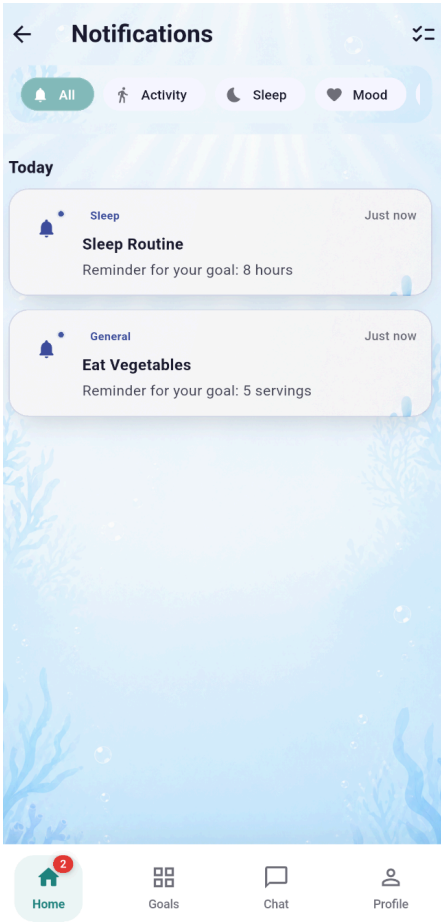

14.1 Notification Categories

Category	What It Includes
All	Complete list of all notifications
Activity	Step goal reminders and activity alerts
Sleep	Sleep routine reminders (e.g., 'Sleep Routine: Reminder for your goal: 8 hours')
Mood	Check-in prompts to log your mood

14.2 Managing Notifications

Each notification card shows:

- Category label (coloured — e.g., blue for Sleep, purple for General)
- Title (e.g. 'Sleep Routine')
- Subtitle (e.g. 'Reminder for your goal: 8 hours')


- Time since received (e.g., 'Just now')



Tap the filter icon (✓≡) in the top-right of the Notifications screen to filter, sort, or mark all notifications as read.

Badge count on the Home navigation icon shows the number of unread notifications.

15 Exporting a Report

ScubaMind can generate a PDF wellbeing report summarising your behavioral data over a chosen time period. This report can be saved to your device or shared with a healthcare provider.

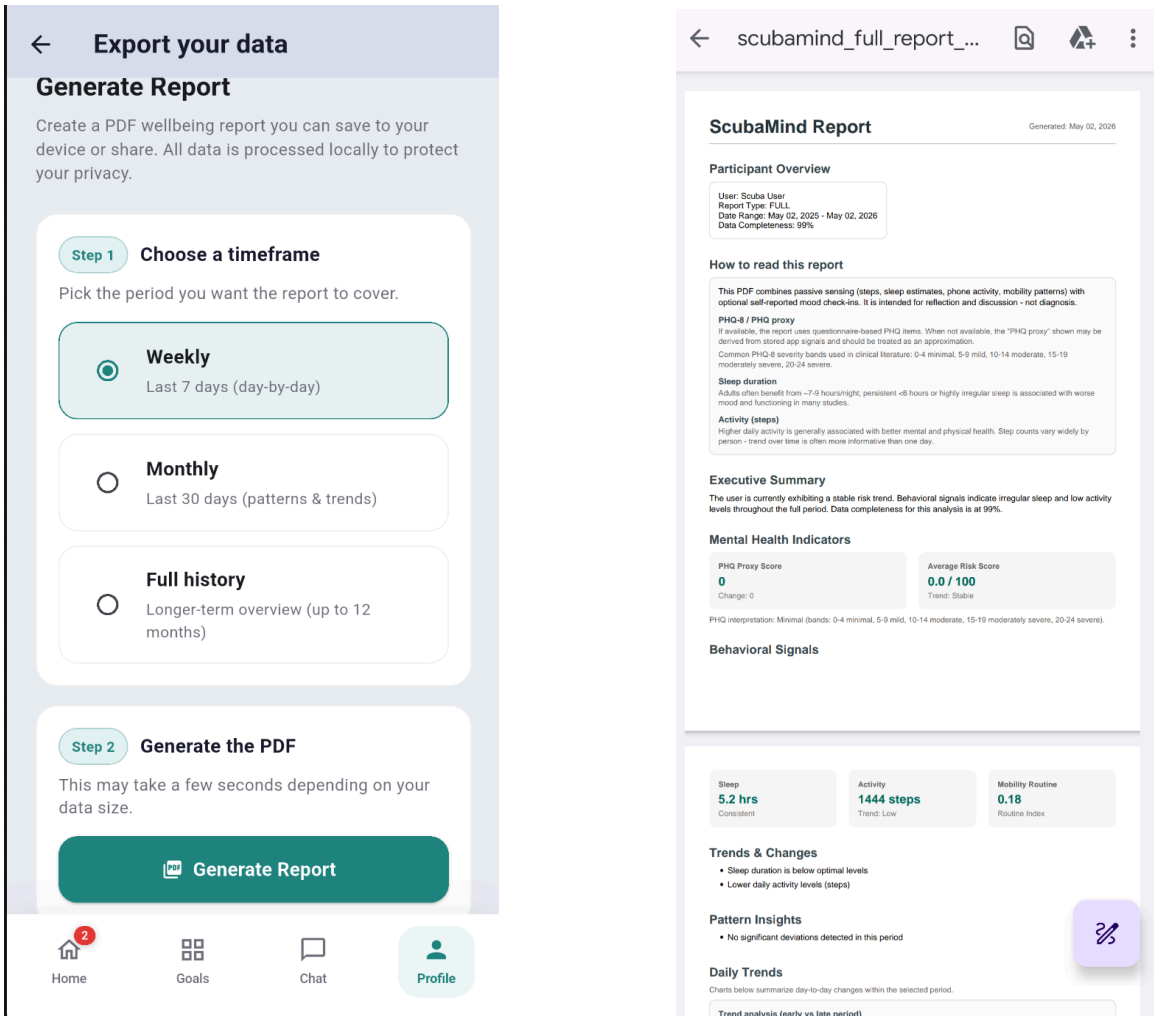

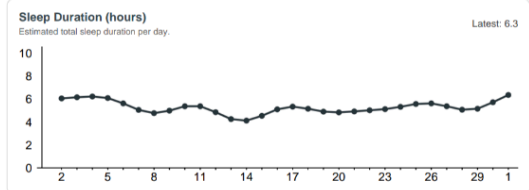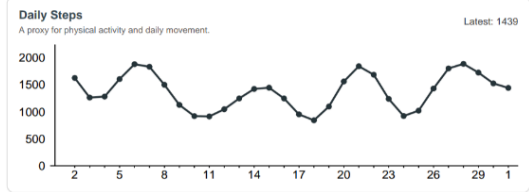

Disclaimer: This report is not a diagnostic tool and does not provide medical advice. It is designed for awareness and support. If your symptoms are severe, persistent, or worsening, consider seeking help from a qualified professional.

Recommendations below are supportive suggestions based on patterns in your data. They are not clinical prescriptions.

- **Sleep:** Your average sleep is below 6 hours. Consider a small, realistic routine change (e.g., a consistent wake time) and track whether mood and energy improve over 1-2 weeks.
- **Activity:** Your steps are on the low side in this window. If it feels doable, try short, light activity (5-10 minutes) on most days and increase slowly.
- **Support idea:** Improve your sleep
- **Support idea:** Increase your daily movement

Tip: If you share this PDF with a clinician, include your recent context (sleep, stressors, schedule changes) for best interpretation.

PHQ-8 severity cutoffs are widely reported in clinical screening literature; sleep duration guidance commonly cites ~7-9 hours/night for adults; physical activity guidelines emphasize regular moderate activity and gradual increases when starting from low levels.

31. Go to Profile → Privacy & Data → Export your data.
32. On the 'Export your data — Generate Report' screen, choose a timeframe:

Option	Description
Weekly	Last 7 days — day-by-day detail
Monthly	Last 30 days — patterns and trends
Full history	Up to 12 months of longitudinal data

33. Tap **Generate Report**. Processing may take a few seconds, depending on your data volume.

34. The PDF opens directly in a system PDF viewer (e.g, 'scubamind_full_report_...pdf').
35. Use the share or save icon in the PDF viewer to save or share the file.

15.2 What the Report Contains

Section	Content
Participant Overview	Username, report type, date range, data completeness percentage
How to Read This Report	Explanations of PHQ score interpretation, sleep duration guidance and , activity guidelines
Executive Summary	Narrative summary of the period (e.g., 'The user is exhibiting a stable risk trend...')
Mental Health Indicators	PHQ proxy score, average risk score, trend direction
Behavioral Signals	Sleep average, activity (steps) average, mobility routine index
Trends & Changes	Bulleted list of notable changes in the period
Pattern Insights	AI-generated pattern observations
Daily Trends Charts	Daily steps chart, sleep duration chart, risk score chart
Recommendations	Personalized supportive suggestions based on patterns
References (brief)	Research-grounded context for interpreting the data
Disclaimer	Clear statement that this is not a diagnostic tool

Sharing with Clinicians

If you share this PDF with a clinician, include your recent context (sleep schedule changes, stressors, life events) for best interpretation.

The report is designed as a discussion aid, not a standalone medical document.

All data in the report is generated locally — no cloud processing is required.

16 Data Security & Privacy

ScubaMind was designed from the ground up with Privacy-by-Design principles. This section explains the technical measures in place to protect your data.

16.1 Local-First Architecture

All sensitive data collection, processing, and storage takes place exclusively on your Android device. This means:

- Your raw GPS coordinates, step counts, sleep data, phone usage statistics, and keyboard text never leave your device.

- Voice recordings and audio files are stored only in the app's private directory (inaccessible to other apps).
- All AI inference (risk score computation, ASR transcription) runs on your device using embedded ML models.
- The app functions fully without an internet connection (except for the AI Chat feature).

16.2 AES-256 Encrypted Database

All behavioral data is stored in an AES-256-encrypted SQLite database (SQLCipher). The database encryption key is:

- Dynamically generated as a UUID at first launch
- Stored in EncryptedSharedPreferences
- Protected by the Android Keystore System (hardware-backed TEE / StrongBox)
- Inaccessible to other apps and to rooted operating systems

16.3 Android Application Sandbox

Android's application sandbox ensures that:

- No other app can access ScubaMind's internal storage, database, voice recordings, or keyboard data without explicit permission.
- Voice recordings are stored in the private app directory with no world-readable permissions
- The Room database files are not accessible to third-party apps.

16.4 Privacy and Security Design Considerations

Standard / Regulation	How ScubaMind Addresses It
KVKK (Turkish Data Protection Law)	Local-first design, explicit consent flows, Data Purge / Right to Erasure feature
GDPR Principles	Data minimization, purpose limitation, user control, right to erasure
OWASP Mobile Security	AES-256 encryption, Keystore key management, private storage for sensitive files
Android Accessibility	Minimum 48×48 dp touch targets, adequate color contrast, screen reader compatibility

17 Troubleshooting & FAQ

Problem	Solution
Home screen shows 'Syncing' for Activity	Activity data requires the Physical Activity permission. Go to Profile → Privacy & Data and ensure Steps & activity is toggled ON. Also check Android Settings → Apps → ScubaMind → Permissions.

Problem	Solution
Sleep data shows 'Quiet for now.'	Sleep tracking depends on screen-on/off events. Ensure you are using your phone normally. For wearable sleep data, ensure Samsung Health and Health Connect are both set up, and Health Connect permissions have been granted.
Voice transcription shows 'Transcribing' indefinitely	The on-device Whisper model can take up to 10 seconds per 30-second clip. If it does not complete after 30 seconds, force-close and relaunch the app. Your recording metadata is preserved in the database.
Risk score shows 0.0 / 100	The system requires at least 1–2 days of collected data to compute a meaningful risk score. Continue using the app, and the score will update as data accumulates.
Chat (Scuby) is not responding.	The AI Chat requires an active internet connection (for the Gemini API). Check your network connection. If connected, the Gemini API may be temporarily unavailable — try again in a few minutes.
App crashes during voice recording	Ensure the Microphone permission is granted. If the crash persists, restart your device and try again. The crash recovery mechanism preserves existing recording metadata.
GPS data showing very limited mobility	Ensure Location permission is set to 'Allow all the time' (not just 'While using the app'). Battery optimization restrictions on some devices may prevent background GPS collection.
Badges are not appearing after completing a goal.	Close and reopen the Goals screen. Badge animations trigger when the goal is marked complete. If the badge does not appear, force-close the app and reopen it.
Data loss after phone restart	All background services are designed to restart automatically after device reboot via the BootReceiver component. If services did not restart, go to Privacy & Data and toggle your preferred signals off and back on.
PHQ questionnaire data not saving	Ensure you tap Submit / Next all the way to the final step. Partially completed questionnaires are not saved. If the problem persists, use Report Bug from the Profile screen.

18 Disclaimers & Important Notices

NOT A MEDICAL DEVICE

ScubaMind is a research and awareness prototype developed as a university Senior Design Project.

It has NOT been validated as a clinical diagnostic tool and has NOT received any medical device certification.

The risk scores, mood indicators, PHQ proxy values and all other outputs are for self-awareness and behavioral reflection only.

They do not constitute a clinical diagnosis, medical advice, or a prescription.

⚠️ AI Model Limitations

AI models trained on population-level data may not generalize equally to all individuals. The models may produce inaccurate or misleading outputs for individuals with behavioral patterns not well-represented in the training data. Models have been trained on limited datasets and require further real-world calibration. Do not make medical decisions based solely on ScubaMind outputs.

Use ScubaMind Responsibly

ScubaMind is intended as a complement to — not a replacement for — professional mental health care.

If you notice concerning trends in your risk score or feel that your mental health is deteriorating, please consult a qualified mental health professional.

Regular PHQ-8 check-ins within the app can help you track changes over time, but a significant increase in your PHQ score should prompt a conversation with a healthcare provider.

The application is designed to encourage seeking professional support, not to discourage it.

19 Contact & Feedback

For technical support, feedback, or bug reports, use the in-app options in Profile → Support, or contact the development team directly.

Contact Method	Details
In-app Feedback	Profile → Support → Feedback
In-app Bug Report	Profile → Support → Report Bug
Email	Profile → Support → Send Email (response within 24 hours)

Team Member	Email
Merve Güleç	merve.gulec@ug.bilkent.edu.tr
Metin Çalışkan	metin.caliskan@ug.bilkent.edu.tr
Muhammed Fatih Başal	fatih.basal@ug.bilkent.edu.tr
Murathan Işık	murathan.isik@ug.bilkent.edu.tr
Yiğit Koşum	yigit.kosum@ug.bilkent.edu.tr

ScubaMind — User Manual

Version 0.1.0 • Bilkent University Senior Design Project — Team T2507 • May 2026

All data processed locally. Your privacy is our priority.